

SpeedyAudit

Quick and Easy Auditing



SpeedyAudit is an easy to use hand hygiene auditing and reporting tool that:

- Automates data collection
- Eliminates transcription
- Provides on-the-spot feedback and reports

You can:

- Capture In/Out, 4 or 5 moments of hand hygiene
- Distinguish between opportunities and indications
- Distinguish between Acute and Long Term Care
- Customize locations and professions
- Capture additional precautions and personal protective equipment
- View reports directly on your mobile device
- Receive compliance rates and data via email
- Integrate the data collected into existing reporting framework

SpeedyAudit Pro:

Gives you access to an online reporting web portal that organizes your data for you, allows you to integrate all data from your facility, create reports and conduct data analysis.



How do I start using SpeedyAudit?





What do I need to get started?

- 1) A device
Apple (iPad, iPhone, iPod) w/ iOS 6+
OR
Android (phone, tablet) w/ JellyBean 4.2.2+
- 2) An email address

Got it, now how do I start auditing?

1. Go to this link in your mobile device browser: <https://speedyaudit.com/speedyaudit>
2. Save the link to your Home Screen
3. Close your browser, and open the link on your home screen.
4. Sign up and Start auditing!

Coming soon to the App Store and Google Play!

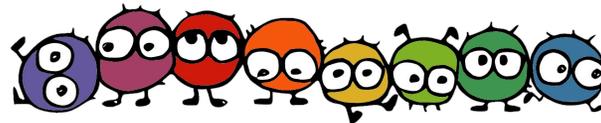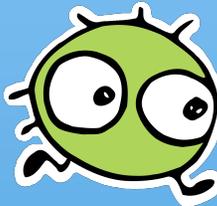

Buttons and Abbreviations



User profile menu: update your email, invite a friend, sign up for SpeedyAudit Pro, and reset your account.



Returns to the main SpeedyAudit Menu.



Cancels the current auditing session.



Ends the current session.



Settings Menu: select your hand hygiene moments, enable/disable track individuals; select facility type, choose your desired interface, choose your language.



Adds a profession / Adds a location.



Functions as 'Next' and saves the data in current session.



Additional Precautions (AP)



Personal protective equipment (PPE)



Tags: include nails, rings, and bracelets



Comments



Reports: view reports



Backs up your data



Help and Support



SpeedyAudit Lite comes with built in tutorials as well as well as a support and Frequently Asked Questions Page.

We welcome feedback, and encourage you to use the Feedback/Comments button in your app to let us know how you like the tool and any areas of improvement you would like to see.

SpeedyAudit Pro comes with more fully featured, personal support.



HANDYMETRICS CORPORATION

203A – 370 QUEENS QUAY WEST, TORONTO, ON M5V 3J3

(917) 426-0586 (TEL/FAX) • SPEEDYSALES@HANDYMETRICS.COM

WWW.HANDYMETRICS.COM



Auditing

SpeedyAudit

- New Session **1**
- Email my Data **2**
- Report **3**
- Info & Help

Tap the "New Session" button to start. SpeedyAudit has 2 auditing interfaces. Select the interface in the menu.

- To audit one healthcare provider at a time use the **Single** interface.
- To audit multiple healthcare providers at the same time, use the **Summary** interface.

1

Summary Interface

Session Summary - Emergency Room

tap here to set profession

Cancel Professions

Physician

Nurse

Medical Student

Social Worker

Physiotherapist

Session Duration: 04:01 0 Opportunities (0 incomplete) with 0 Indications

After choosing the Auditing Location, select the HCP profession by tapping on "Add HCP", and then tapping on "tap here to set profession".

You may scroll down the list to find more professions or tap on to add a profession. When you are done selecting a profession, tap on the button.

Session Summary - Emergency Room

Physician

BEF-P/E BEF-ASP AFT-BFL AFT-P/E

Missed Wash Rub

16s Missed Wash Rub

16s Missed Wash Rub

BEF-P/E BEF-ASP AFT-BFL AFT-P/E

Missed Wash Rub

tap here to set profession

Cancel Professions

Physician

Nurse

Medical Student

Social Worker

Physiotherapist

Session Duration: 01:01 3 Opportunities (0 incomplete) with 4 Indications

Each line represents an opportunity. Select the moment(s) and indicate whether these were missed or complied. When tapping on Rub or Wash a timer will appear.

You may add more opportunities to the existing HCP by tapping on . Here you may indicate if the patient was under Additional Precautions (AP), the PPE used, and additional information (e.g. whether the HCP is wearing rings or bracelets), and you may enter any other comments relevant to the hand hygiene opportunity.

To add another profession to the current session, tap on "+Add HCP", and then "tap here to set profession". To end the session tap on .

1

Single interface

Professions

Location: Emergency Room

Patient Transporter

Radiology Tech

Respiratory Therapist

Dietician

PSA, PSW, PCA

Other

+ Add Custom Profession

After choosing the Auditing Location, select the HCP profession by tapping on an item from the list. To view the rest of the menu, drag the list up or down.

Opportunity

Location: Emergency Room
Profession: Radiology Tech

Additional Precautions

Contact Droplet Airborne

PPE

Gloves Gown Goggles

Mask N95

Compliance Precautions Tags

To select Additional Precaution (AP) or if the HCP is donning any Personal Protective Equipment (PPE), open the page by tapping on 'Precautions' located on the tab bar, below.

To select an AP or PPE, tap on the appropriate button. Tap on the button again to deselect the item. You can select multiple APs / PPEs.

Opportunity

Location: Emergency Room
Profession: Radiology Tech

Missed Wash Rub

Session Duration: 00:17 0 Opp - 0 Ind

Compliance Precautions Tags

To select an indication, tap on an arrow. Tap the arrow again to deselect the item. When appropriate, you can select multiple indications.

Tap 'Missed,' 'Wash,' or 'Rub' to indicate the compliance of the current opportunity.

Opportunity

Location: Emergency Room
Profession: Radiology Tech

Tags

Nails Rings Bracelets

Comments

Compliance Precautions Tags

To indicate if the patient has inappropriate nails, rings or bracelets, open the page by tapping on 'tags' located on the tab bar.

You may also add relevant comments about the audit by tapping on the white text box under 'Comments.'

To finish recording an opportunity and start a new one, tap on .

2

Emailing your data

Email my Data

Sep 22, 2014 - Emergency Room
1 Opportunity - 1 Indication Send

Sep 22, 2014 - Emergency Room
1 Opportunity - 2 Indications Send

Sep 18, 2014 - Intensive Care Unit
3 Opportunities - 7 Indications Send

Once you've ended your session, you can send your data via email. To send your data, tap on 'Email my Data' from the SpeedyAudit main menu. Before sending your data you will need to verify your email address. Once you have verified your email address, tap 'Send.' You will then receive an email with the observations you collected in a CSV format.

3

Creating a Report



To View a report, click on the 'Report' button in your main menu. The report will show aggregate data from all of your auditing sessions.

If you view a report during an auditing session, then the report will only show data from the current auditing session.